



Opinion

Volume 7 Issue 2 – september 2017  
DOI: 10.19080/ARGH.2017.07.555709

Adv Res Gastroentero Hepatol

Copyright © All rights are reserved by Hyunghun Kim

# Side Attachment Method: Choice for Easy and Safe Intra-gastric Balloon Insertion



Hyunghun Kim\*

Department of Gastroenterology, Kangnam Peter's Hospital, South Korea

Submission: September 11, 2017; Published: September 25, 2017

\*Corresponding author: Hyunghun Kim, Department of Gastroenterology, Kangnam Peter's Hospital, Seoul, Korea Address: Kangnam Peter's hospital, 2633 Nambusunhwan-ro, Gangnam-gu, Seoul, South Korea, Tel: +82-2-1544-7522; Fax: +82-2574-9414; Email: relyonlord@gmail.com

## Abstract

The intra-gastric balloon is an effective method for weight reduction in obese people. The common problem of intra-gastric balloon insertion is a blind performance. Practitioners cannot see the insertion from oral cavity to the stomach clearly. To overcome this point, side attachment method was developed. The introduction kit is attached to the shaft of the endoscope side by side, so the introduction kit does not cover the camera of an endoscope. Side attachment method provides clear field of vision during insertion from oral cavity to the stomach, same to diagnostic endoscopy.

**Keywords:** Intra-gastric balloon; Endoscopy; Obesity

**Abbreviation:** IGB, intra-gastric balloon; SAM, side attachment method

## Introduction

The intra-gastric balloon (IGB) technique has become an effective and safe method of achieving significant weight reduction in obese people. IGBs allow patients to feel fullness and finally decrease their food intake. IGBs replace gastric luminal volume and may distend the stomach, potentially inducing neurohormonal effects and changes in motility [1-5]. Most IGBs are blindly inserted into the stomach without endoscopic guide, finger technique [6]. End-Ball® (Endalis, Brignais, France) uses endoscope as a guide for passage of the balloon into the stomach; however, the endoscopic view is nearly blind because an introduction kit covers the camera of an endoscope, front attachment method. The common problem of these procedures is a blind performance. If a practitioner can insert an introduction kit with diagnostic visual field, insertion procedure will be easier and safer; blind insertion can damage oral cavity and gastrointestinal tract involuntarily.

## Side Attachment Method

Side Attachment Method (SAM) was developed to provide clear visual field during insertion. In SAM, the introduction kit is attached to the shaft of the endoscope side by side, so the introduction kit does not cover the camera of an endoscope. Additionally, SAM with J turn allows a practitioner to notice that an introduction kit is in the stomach. SAM with J turn also enables a practitioner to watch inflation process time to time [8].

## Practice of Side Attachment Method

- i. A 3.5cm X1cm sized window was made at the distal part of the pouch for releasing a balloon before side attachment to an endoscope.



**Figure 1:** An introduction kit is attached to the endoscope at 2 cm, 20 cm, 30 cm from the end of the endoscope. The window side is marked with dotted line.

- ii. The endoscope was then connected to the introduction kit by side attachment method (Figure 1). The opposite side of the window was attached to the scope by 2.5cm width paper plaster, so the window is placed at the opposite direction of attached side. The tip of the pouch was attached to the point 2cm distal to the end of the scope, and the catheter was attached to the points 20cm and 30cm distal to the end of the endoscope.

iii. Inserting the endoscope with introduction kit as a diagnostic endoscopy insertion. While passing through the oral cavity, the pyriform sinus, the esophagus, all structures were clearly observed as diagnostic endoscopy.

iv. After inserting an introduction kit into the stomach, J turn was used to confirm whether the introduction kit is in the stomach (Figure 2a).

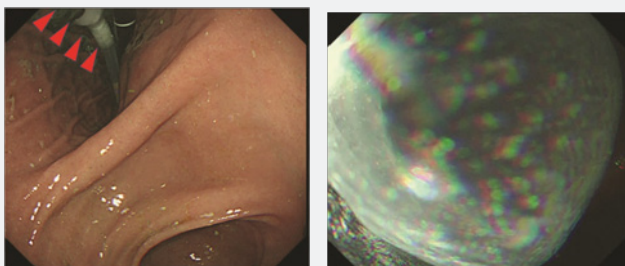


Figure 2a


Figure 2b

**Figure 2:** J turn of an endoscope during side attachment method implantation (a) J turn of the endoscope reveals that the pouch is placed in the stomach (red arrow heads). (b) J turn of the endoscope shows an inflating balloon.

v. The balloon began to be filled. When 100 cc air and 150 cc saline were applied, the balloon was pushed out of the pouch; J turn also showed this process (Figure 2b) and filled with saline to its designated volume (600 cc) under endoscopic control.

vi. After then, a feeding catheter was pulled up with an endoscope at the cardia, and the IGB was disconnected from the feeding catheter.

vii. Finishing implantation, follow-up endoscopy revealed that the IGB was well placed in the fundus (Figure 3).

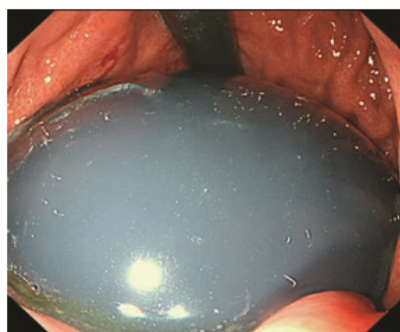

**Figure 3:** An intragastric balloon is well located in the stomach.

### Conclusion

SAM is a very useful technique for implanting IGB. Ordinary front visual field during insertion is a pivotal advantage of SAM. Conventional technique caused nearly blind insertion, so practitioners cannot but rely on tactile sensation instead of visual sensation. SAM is a reasonable choice for easy and safe IGB implantation.

### References

1. Mathus-Vliegen EM (2014) Endoscopic treatment: the past, the present and the future. *Best Pract Res Clin Gastroenterol* 28(4): 685-702.
2. Geliebter A, Westreich S, Gage D (1988) Gastric distention by balloon and test-meal intake in obese and lean subjects. *Am J Clin Nutr* 48(3): 592-594.
3. Geliebter A, Melton PM, McCray RS, Gage D, Heymsfield SB, et al. (1991) Clinical trial of silicone-rubber gastric balloon to treat obesity. *Int J Obes* 15(4): 259-266.
4. Mathus-Vliegen EM, Eichenberger RI (2014) Fasting and meal-suppressed ghrelin levels before and after intragastric balloons and balloon-induced weight loss. *Obes Surg* 24(1): 85-94.
5. Mathus-Vliegen EM, de Groot GH (2013) Fasting and meal-induced CCK and PP secretion following intragastric balloon treatment for obesity. *Obes Surg* 23(5): 622-633.
6. Stimac D, Majanovic SK, Turk T, Kezele B, Licul V, et al. (2011) Intragastric balloon treatment for obesity: results of a large single center prospective study. *Obes Surg* 21(5): 551-555.
7. Buzga M, Kupka T, Siroky M, Narwan H, Machytka E, et al. (2016) Short-term outcomes of the new intragastric balloon End-Ball(R) for treatment of obesity. *Wideochir Inne Tech Maloinwazyjne* 11(4): 229-235.
8. Hyunghun Kim (2017) <https://youtu.be/DaCEXaz1Ofk>



This work is licensed under Creative Commons Attribution 4.0 License  
DOI: [10.19080/ARGH.2017.07.555709](https://doi.org/10.19080/ARGH.2017.07.555709)

### Your next submission with JuniperPublishers will reach you the below assets

- Quality Editorial service
- Swift Peer Review
- Reprints availability
- E-prints Service
- Manuscript Podcast for convenient understanding
- Global attainment for your research
- Manuscript accessibility in different formats  
( Pdf, E-pub, Full Text, audio)
- Unceasing customer service

Track the below URL for one-step submission  
<https://juniperpublishers.com/online-submission.php>